




Coag **L**

THE COMPACT AUTOMATED COAGULOMETER FOR MEDIUM SIZE LABORATORIES

4 channels, 160 PT test/hour*
continuous random access
coagulation, chromogenic and immunology tests

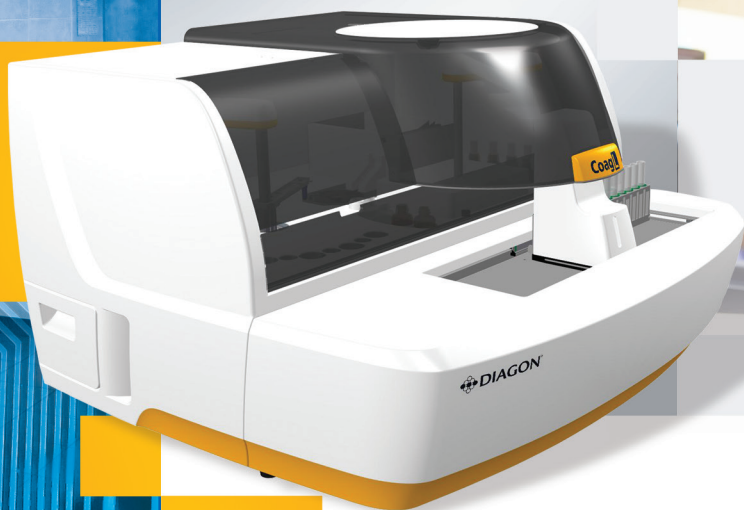

CoagL

AUTOMATED COAGULATION INSTRUMENT

Key features:

- Random access of optical clotting, immunologic and chromogenic assays
- 4 independent dual-wavelength measuring channels
- Increased throughput: up to 160 PT/hour*
- Master curve calibration for routine tests and D-dimer, auto dilution, automatic reflex / rerun feature
- Filtering algorithm to avoid lipid interference
- Tilted reagent holder for reduced dead volumes
- Multilot use of reagents
- Integrated barcode reader
- Optional 24/7 real-time remote access
- STAT sample positions
- Walkaway capacity: 70 samples / 500 cuvettes

Hardware:

- External colour LCD touchscreen
- 14 cooled + 4 reagent positions
- Unique cuvette transport system
- Bidirectional LIS connection

Software:

- Innovative user interface
- Multilanguage software
- Smart QC menu
- Extended user access level and password protection

*measured using Diagon Dia-PT Recombinant test

TECHNICAL PARAMETERS:

Dimensions:	Height 510 mm; Width 650 mm; Depth 750 mm
Weight:	80 kg
Power supply:	AC 100-250 V
Data output:	6 USB interface, network interface
Incubation block:	37 °C ± 0,1 °C
Wavelength detection:	405 nm, 570 nm and 640 nm



EUROPE

Diagon Ltd - Hungary
1047 Budapest, Baross u. 48-52.
phone.: (+36 1) 369 6500
fax: (+36 1) 369 6301
e-mail: diagon@diagon.com

Diagon GmbH – Austria
GZO – Zukunftsweg 1,
Siegendorf, Burgenland 7011
phone: (+43 2687) 42090
e-mail: diagonmbh@diagon.com

www.diagon.com

SOUTH AMERICA

Diagon Ltda - Brasil
Av. Joaquim José Diniz, 787
CEP 31910-520
Bairro Fernão Dias
Belo Horizonte/MG
phone: (+55 31) 3488 6777
fax: (+55 31) 3488 6778
e-mail: diagon@diagon.com.br
www.diagon.com.br

ASIA

PT Inti Diagonama Selaras - Indonesia
Perkantoran Plaza Pasifik A3/59 -
Jl. Raya Boulevard Barat
Kelapa Gading Permal Jakarta 14250
phone.: (+62 21) 458 40649, 458 40650
fax: (+62 21) 458 45024
e-mail: dsu@centrin.net.id